

## Drugs of Abuse Testing



# Target Screening and Quantitative Confirmation

## **MassTox<sup>®</sup> Drugs of Abuse Testing in Urine by LC-MS/MS**

Target screening and quantitative confirmation of more than  
100 drugs and metabolites within the same run

# LC-MS/MS Urine Assay for Both: Target Screening and Quantitative Confirmation

- Assay includes 108 drugs and metabolites
- Target screening and quantitative confirmation within the same run
- Run time of 12 minutes
- Simple sample prep combined with high selectivity and sensitivity
- Internal standards for 98 analytes

The LC-MS/MS based assay **MassTox**<sup>®</sup> Drugs of Abuse Testing offers the target screening of more than 100 drugs and metabolites with an instrument-dependent analysis time of approximately 12 minutes. A quantitative confirmation can also be performed within the same run. The CE-IVD validated solution includes everything, eliminating the need for developing and maintaining several in-house methods or for preparing QC materials for a multitude of analyte compounds.


It also provides a significantly higher selectivity and sensitivity than current immunoassays and is designed to work below cut-offs defined by GTFCh and EWDTs. The assay can be used for target screening, quantitative confirmation, or both. 98 analytes are safeguarded by isotopically labelled internal standards. The high stability of the **6PLUS1**<sup>®</sup> multilevel urine calibrator set and three **MassCheck**<sup>®</sup> urine controls ensure reliable results with highest precision.

## Parameter Menu

|   |  |  |  |
|---|--|--|--|
| <b>Amphetamines</b><br>Amphetamine<br>BDB<br>Butylone<br>2C-B<br>2C-I<br>Cathinone<br>MBDB<br>MDA<br>MDEA<br>MDMA<br>MDPV<br>Mephedrone<br>Methamphetamine<br>Methaqualone<br>Methylone<br>Methylphenidate<br>PMA<br>Ritalinic Acid | <b>Benzodiazepines</b><br>Alprazolam<br>7-Aminoclonazepam<br>7-Aminoflunitrazepam<br>7-Aminonitrazepam<br>Bromazepam<br>Brotizolam<br>Chlordiazepoxide<br>Clobazam<br>Clonazepam<br>Demoxepam<br>Desalkylflurazepam<br>Desmethylflunitrazepam<br>Diazepam<br>Estazolam<br>Flunitrazepam<br>Flurazepam<br>Lorazepam<br>Lormetazepam<br>Medazepam<br>Midazolam<br>Nitrazepam<br>Norclobazam<br>Nordiazepam<br>$\alpha$ -OH-Alprazolam<br>3-OH-Bromazepam<br>$\alpha$ -OH-Midazolam<br>$\alpha$ -OH-Triazolam<br>Oxazepam<br>Prazepam<br>Temazepam<br>Triazolam | <b>Booster</b><br>Gabapentin<br>Pregabalin<br>Promethazine<br>Quetiapine       | <b>Opiates/Opioids</b><br>Acetylcodeine<br>Buprenorphine<br>Codeine<br>Dihydrocodeine<br>EDDP<br>Fentanyl<br>Hydrocodone<br>Hydromorphone<br>Meconin<br>Meperidine<br>Methadone<br>6-Monoacetylmorphine<br>Morphine<br>Naloxone<br>Naltrexone<br>Norbuprenorphine<br>Norcodeine<br>Norfentanyl<br>Normeperidine<br>Nortapentadol<br>Nortilidine<br>O-Desmethyltramadol<br>Oxycodone<br>Oxymorphone<br>Papaverine<br>Propoxyphene<br>Sufentanil<br>Tapentadol<br>Thebaine<br>Tilidine<br>Tramadol |
| <b>Barbiturates</b><br>Allobarbitol<br>Amobarbitol<br>Barbitol<br>Butalbital<br>Hexobarbitol<br>Pentobarbitol<br>Phenobarbitol<br>Secbutabarbitol<br>Secobarbitol<br>Thiopental   |  | <b>Cannabinoids</b><br>THC-COOH  |  |
|   |  | <b>Cocaine</b><br>Benzoylcegonine<br>Cocaethylene<br>Cocaine<br>Norcocaine     |  |
|   |  | <b>Others</b><br>Ketamine<br>LSD<br>Mescaline<br>Norketamine<br>O-H-LSD<br>PCP |  |
|   |  | <b>Z-Drugs</b><br>Zaleplon<br>Zolpidem<br>Zopiclone                            |  |

# Sample Preparation: Fast and Easy

## Sample Preparation in Reaction Vials\*

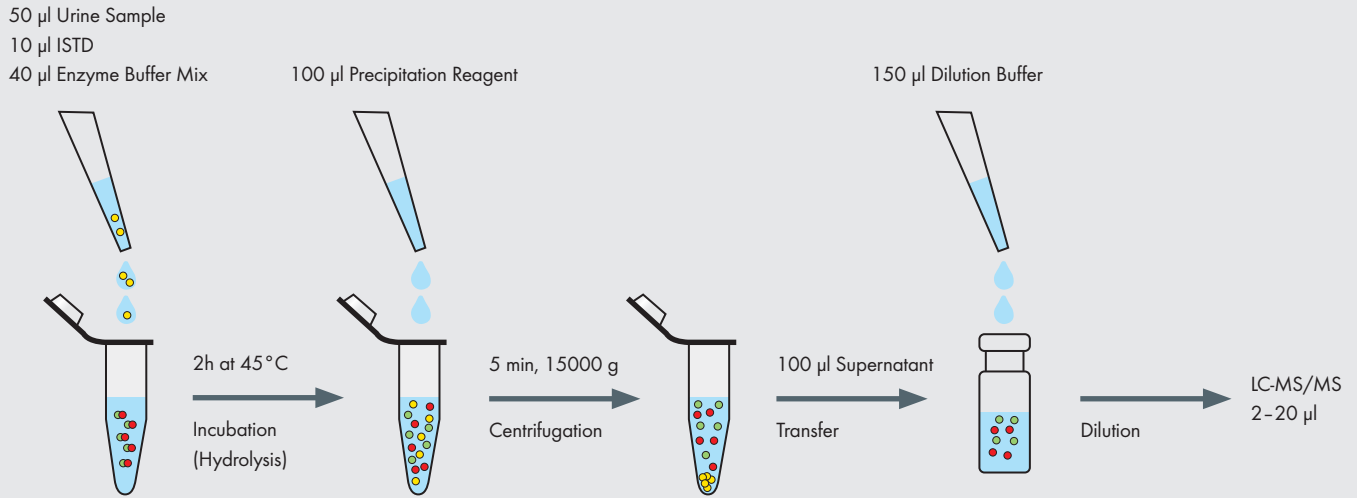

### Pre-analytic Treatment

|           |   |
|-----------|---|
| Specimen  | urine   |
| Stability | according to the Swiss Guidelines Committee for Drugs of Abuse Testing (SCDAT), samples are stable for seven days at +2 to +8 °C and for six months at -18 °C |

### Specifications


|                            |                                    |
|----------------------------|------------------------------------|
| Linearity for all analytes | from 0.25 to 30 x of cut-off value |
| Limit of quantification    | down to 13 ng/l                    |
| Intraassay                 | CV = 1-9 %                         |
| Interassay                 | CV = 1-11 %                        |
| Recovery                   | 85-115 %                           |
| Analysis time              | 12 min                             |

\*Assay also available for 96 Well Plates

## 100% Hydrolysis – 0% Doubt

Many drugs undergo glucuronidation prior to excretion, which requires an enzymatic or acidic hydrolysis – a challenge for most assays and in-house methods<sup>[1]</sup>. In contrast, the **MassTox**® Drugs of Abuse Testing assay from Chromsystems includes a carefully selected


enzyme that has proven to function with all 108 drugs, including codeine and other substances usually difficult to hydrolyse. This ensures a complete and selective hydrolysis of all drugs covered in the assay.

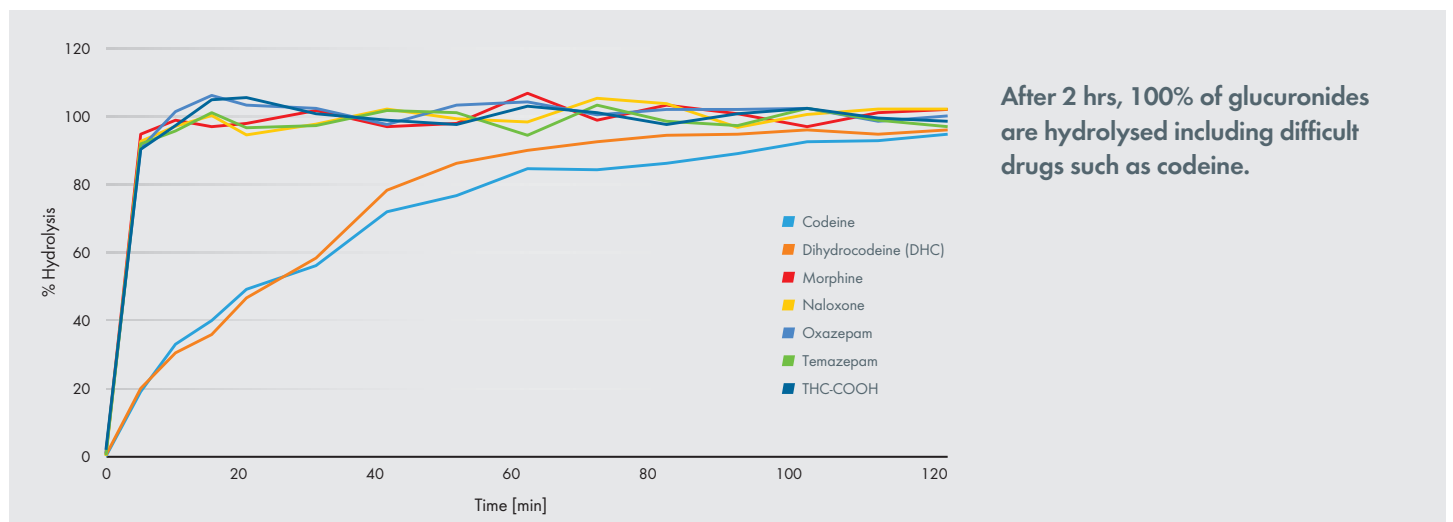

## Proven Accuracy and Reliability

A lab experienced in routine drugs of abuse testing performed a comparative analysis of the **MassTox**® assay with an LC-MS/MS based in-house method. Values obtained correspond very well with each other at concentrations within the linearity of the method. This demonstrates

the accuracy of data and the suitability of the Chromsystems assay to replace LC-MS/MS based in-house methods in the forensic and clinical laboratory. The **MassTox**® assay also has been successfully used in proficiency testing schemes of GTFCh, Instand and RfB.

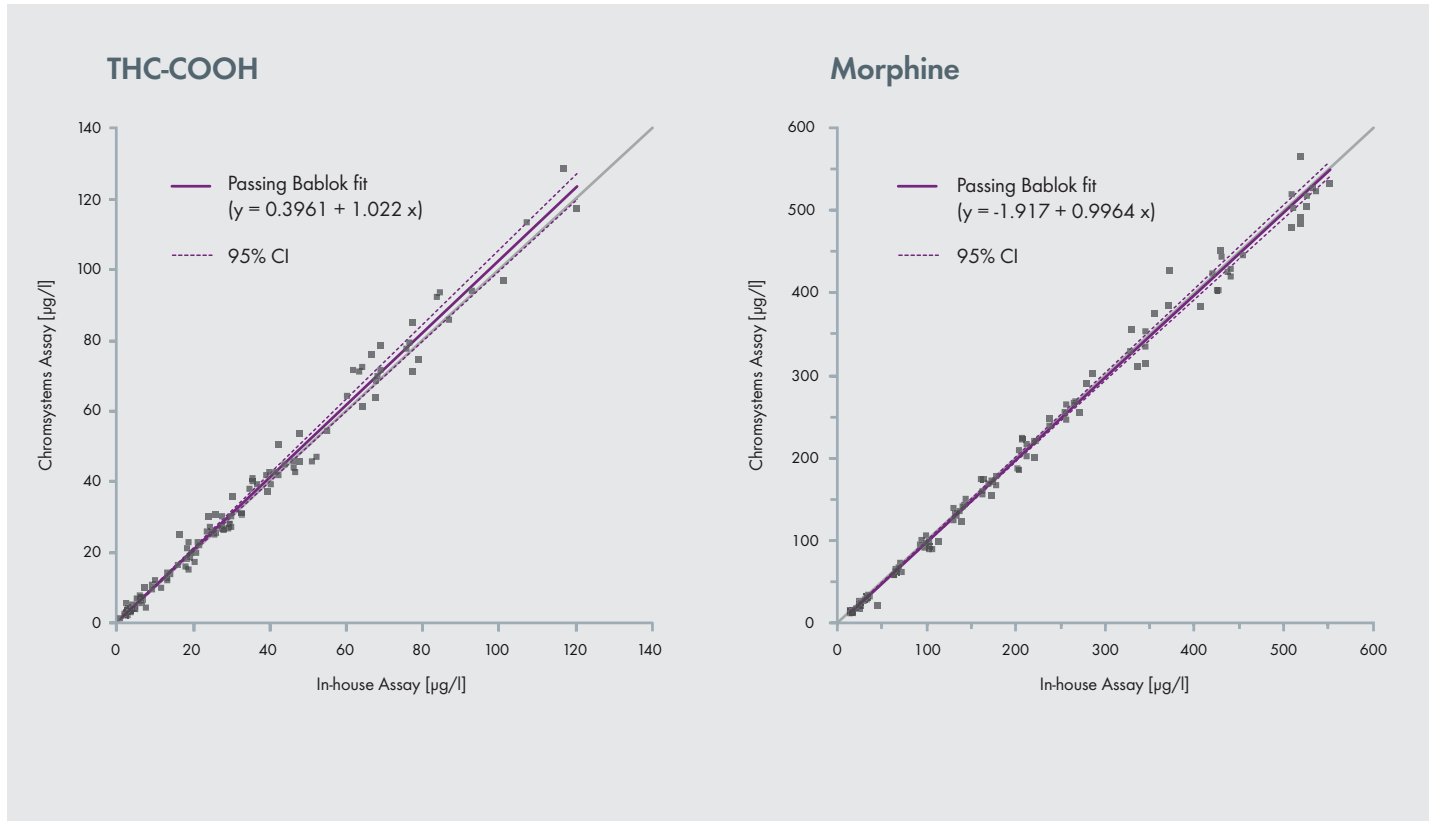


## LC-MS/MS: Faster Sample Prep – Lower Workload

A comparison of the utility of the LC-MS/MS assay in the daily routine with GC-MS also has been made. Results not surprisingly demonstrate that LC-MS/MS allows the direct analysis of the specimens and eliminates the need for derivatisation steps. This makes the sample prep much simpler and less time-consuming in compar-

ison to GC-MS. Another advantage is the need for much lower sample volumes. This means: when comparing our LC-MS/MS assay with GC-MS, then it provides comparable data, but with the advantage for LC-MS/MS to offer a faster sample prep and lower sample volumes.

|   | LC-MS/MS    | GC-MS        |
|---|-------------|--------------|
| Total run time                          | 12 min      | 10–15 min    |
| Sample prep time (excluding hydrolysis) | 2 h per day | 6 h per day  |
| Sample volume                           | 50 µl       | 1000–5000 µl |


## 108 Drugs/Metabolites – 1 Analytical Run

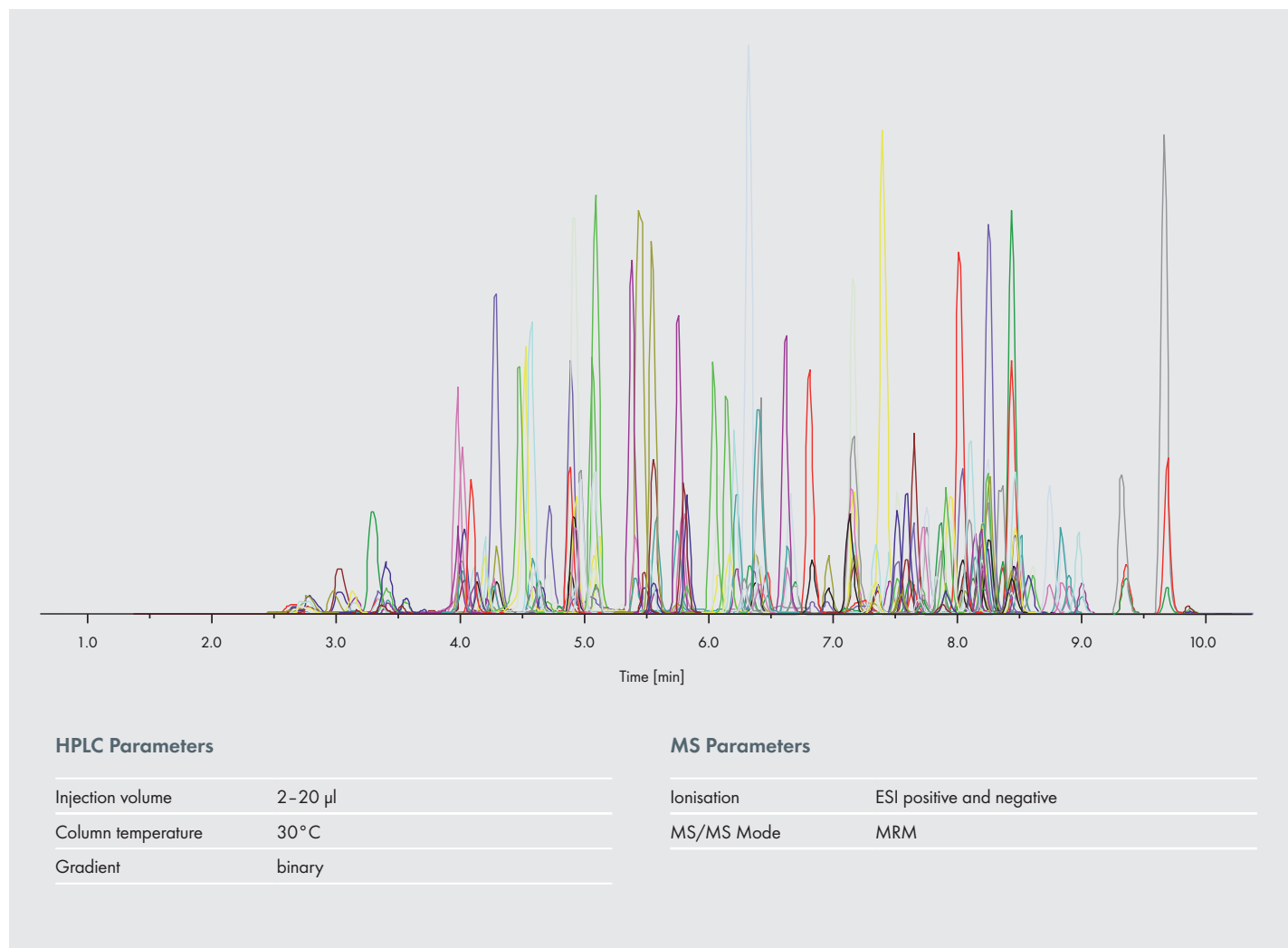

## Calibrator and Control Material

The Chromsystems assay is provided with a 6PLUS1® Multilevel Urine Calibrator Set and **MassCheck**® controls with three levels that cover the relevant concentration range, including cut-off values from GTFCh and EWDTs. The material includes a lyophilised urine

matrix that is reconstituted with a buffer. It has been developed to provide highest stability for all 108 analytes. For target screening alone, a Urine Screening Standard Set is also available.

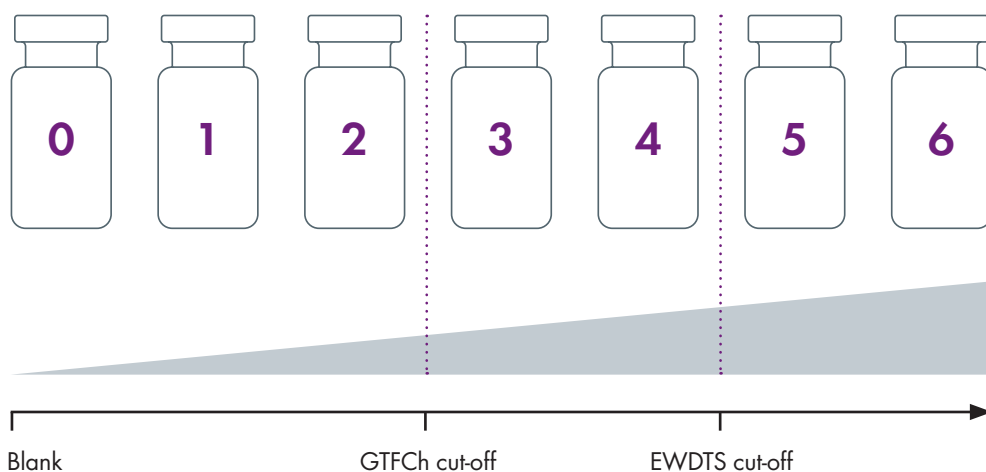

The Chromsystems assay covers the cut-off values defined by GTFCh and EWDTs.

## MassTox<sup>®</sup> Drugs of Abuse Testing

**96000** **MassTox<sup>®</sup> Drugs of Abuse Testing in Urine (with Reaction Vials) – For 400 tests**

**96000/WP** **MassTox<sup>®</sup> Drugs of Abuse Testing in Urine (with 96 Well Plates) – For 400 tests**

### Components available separately

- 96001 Mobile Phase A, 700 ml
- 96002 Mobile Phase B, 700 ml
- 96005 Precipitation Reagent, 40 ml
- 96009 Rinsing Solution, 500 ml
- 96010 Dilution Buffer, 30 ml
- 96046 Internal Standard Set, consisting of:
  - Internal Standard Mix (lyoph.), 4 x 1 ml
  - Reconstitution Buffer, 5 ml
- 96078 Enzyme Solution Set, consisting of:
  - Enzyme Reagent (lyoph.) 2 x 8 ml
  - Hydrolysis Buffer, 20 ml
- 33006 Reaction Vials, 100 pcs. (for 96000)
- 96057 96 Well Plates, 5 pcs. (for 96000/WP)
- 96058 Collection Plates, 5 pcs. (for 96000/WP)
- 96059 Pierceable Adhesive Seals for 96 Well Plates, 5 pcs. (for 96000/WP)

### Startup Accessories

- 96100 Analytical Column, equilibrated, with test chromatogram, 1 pc.
- 96033 Urine Screening Standard Set (lyoph.), 4 x 1 ml
- 96012 System Check Solution, 1 ml
- 96015 Tuning Mix 1, Analytes and Internal Standards, 1 ml
- 96016 Tuning Mix 2, Analytes and Internal Standards, 1 ml
- 96017 Tuning Mix 3, Analytes and Internal Standards, 1 ml
- 15090 UHPLC Stainless Steel Prefilter Housing, 1 pc.
- 15091 Stainless Steel Prefilter, 0.2 µm, 3 pcs.

### Multilevel Calibrator and Controls

- 96040 6PLUS1<sup>®</sup> Multilevel Urine Calibrator Set (lyoph.), 7 x 1 ml
- 0410 **MassCheck<sup>®</sup>** Drugs of Abuse Testing Urine Control Set, Level I (lyoph.), 4 x 1 ml
- 0420 **MassCheck<sup>®</sup>** Drugs of Abuse Testing Urine Control Set, Level II (lyoph.), 4 x 1 ml
- 0430 **MassCheck<sup>®</sup>** Drugs of Abuse Testing Urine Control Set, Level III (lyoph.), 4 x 1 ml

### Hydrolysis Products

- 0470 **MassCheck<sup>®</sup>** Drugs of Abuse Testing Urine Hydrolysis Control Set (lyoph.), 4 x 1 ml
- 96078 Enzyme Solution Set, consisting of:
  - Enzyme Reagent (lyoph.) 2 x 8 ml
  - Hydrolysis Buffer, 20 ml

These products are not available in all countries. Please contact your local Chromsystems sales representative for availability in your region.

Please note: The product information provided shows typical values as an example only. Please review the information provided in instruction manual and product information leaflets.

Reference: [1] Wang et al., Incomplete Recovery of Prescription Opioids in Urine using Enzymatic Hydrolysis of Glucuronide Metabolites. J. Analytical Toxicology, (2019), 571-575.